

63X\_1\_2\_525

### Microscope info:

| Image                            |                      | GreenBeads63x.lif-Series004_bead395 |                         |            |                       |                       |                     |  |
|----------------------------------|----------------------|-------------------------------------|-------------------------|------------|-----------------------|-----------------------|---------------------|--|
| image's date 2024-10-26 09:47:32 |                      |                                     |                         |            |                       |                       |                     |  |
| creation                         | method used          | from file                           | from file creation date |            |                       |                       |                     |  |
| Actual image                     | 12                   |                                     |                         |            |                       |                       |                     |  |
| Microsco                         | ope type             | WideField                           |                         |            |                       |                       |                     |  |
|                                  | NA                   | 1.2                                 | 1.2                     |            |                       |                       |                     |  |
| Objective                        | im. refractive index | 1.333                               |                         |            |                       |                       |                     |  |
|                                  |                      | Wavel                               | engths                  |            | sampling (X,Y,Z)      |                       |                     |  |
| Chan                             | nel(s)               | Ex.<br>(nm)                         | Em.<br>(nm)             | Saturation | Nyquist (µm)          | Found (µm)            | Nyquist/found ratio |  |
| Channel 0                        |                      |                                     | 525.0                   | none       | 0.109x0.109x0.<br>349 | 0.103x0.103x0.<br>099 | 0.9, 0.9, 0.3       |  |
| Bead original pix                | 1071.0,              | 1690.0                              |                         |            |                       |                       |                     |  |

## Warnings:

(No saturated pixels detected). (All channels sampled following Shannon-Nyquist criterion). (A subresolution bead is used for all channels).

#### Resolution table:

| Channel                 | Sig/Backgn<br>d ratio | Dimension | Measured<br>FWHM<br>(µm) | theory (µm) | Fit<br>Goodness | Mes./theory ratio |
|-------------------------|-----------------------|-----------|--------------------------|-------------|-----------------|-------------------|
|                         |                       | Х         | 0.274                    | 0.223       | 1.0             | 1.23              |
| Channel 0 (em. 525.0nm) | 1.3                   | Υ         | 0.28                     | 0.223       | 0.99            | 1.25              |
| 323.01111)              |                       | Z         | 0.891                    | 0.86        | 0.98            | 1.04              |

Green: within specifications, red: outside specifications (ie. XY ratios above 1.5 or Z ratio above 2.0)

#### Lateral asymmetry ratios:

| Channel                 | Ratio |
|-------------------------|-------|
| Channel 0 (em. 525.0nm) | 0.98  |


## Detailed channel detection info:

## Channel #0





| Channel 0 (em. 525.0nm) |      |           |       |              |  |
|-------------------------|------|-----------|-------|--------------|--|
| Sig./Backgnd ratio      | LAR  | Dimension | FWHM  | Fit goodness |  |
|                         |      | X         | 0.274 | 1.0          |  |
| 1.3                     | 0.98 | Υ         | 0.28  | 0.99         |  |
|                         |      | Z         | 0.891 | 0.98         |  |

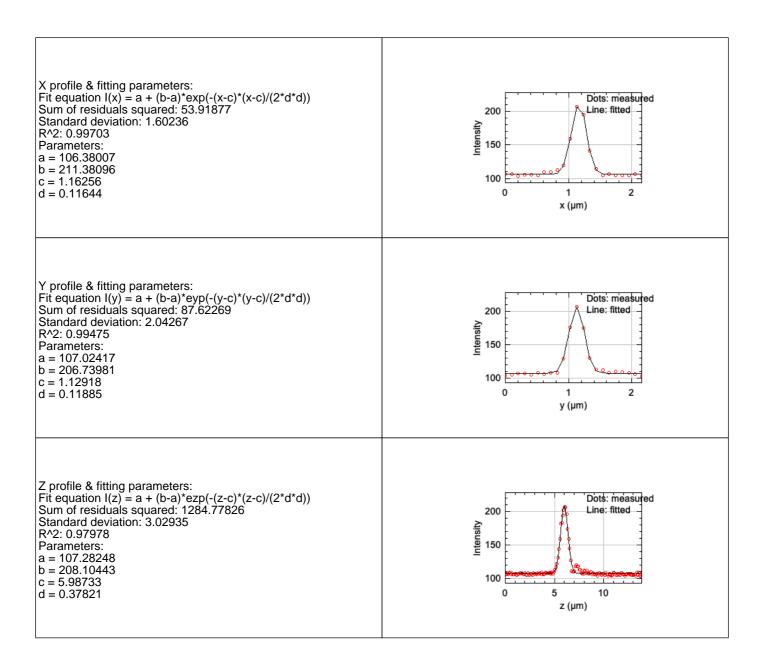

### Analysis parameters

|                                 | Tool                                                   | PSF Profiler (batch)                                                                                                                              |  |  |
|---------------------------------|--------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------|--|--|
| Tool &<br>Operator              | Versions                                               | MetroloJ_QC v1.3.1.1, ImageJ v2.14.0/1.54f, Java v1.8.0_322, OS Mac OS X                                                                          |  |  |
|                                 | Operator & date                                        | SO, October 31, 2024 8:20 AM                                                                                                                      |  |  |
| data                            | result folder                                          | /Users/oggsc/Documents/OM/ImageAnalysis/QC/Thunder/<br>PSF/20241015/63X_1_2_525//Processed/63X_1_2_525/Gr<br>eenBeads63x.lif - Series004/bead395/ |  |  |
| data                            | Type of saved data                                     | .pdf, .jpg, .xls                                                                                                                                  |  |  |
|                                 | Input data bit depth                                   | 12                                                                                                                                                |  |  |
| Dim                             | ension order                                           | XY-(C)Z                                                                                                                                           |  |  |
| Discard s                       | aturated samples                                       | true                                                                                                                                              |  |  |
|                                 | Bead detection threshold                               | Legacy                                                                                                                                            |  |  |
|                                 | Center detection method                                | Legacy Maximum Intensity                                                                                                                          |  |  |
|                                 | Discard bead if more than one particle are thresholded | true                                                                                                                                              |  |  |
|                                 | Background annulus thickness in µm                     | 0.5                                                                                                                                               |  |  |
| Beads                           | Background annulus<br>distance to bead edges<br>in µm  | 0.5                                                                                                                                               |  |  |
| Deads                           | Multiple beads in image                                | true                                                                                                                                              |  |  |
|                                 | Bead identification method                             | Using Find Maxima (prominence of 50.0)                                                                                                            |  |  |
|                                 | Bead size (µm)                                         | 0.1                                                                                                                                               |  |  |
|                                 | Bead crop Factor                                       | 5.0                                                                                                                                               |  |  |
|                                 | Cropped ROI size in µm                                 | 2.31x2.31 (using bead size & background annulus parameters)                                                                                       |  |  |
|                                 | Bead rejection distance to top/bottom                  | 2.0 μm                                                                                                                                            |  |  |
| Square Root PSF Image displayed |                                                        | true                                                                                                                                              |  |  |
| Tolerance                       | Applied in this report                                 | true                                                                                                                                              |  |  |
|                                 | X & Y FWHM ratios valid if below                       | 1.5                                                                                                                                               |  |  |
|                                 | Z FWHM ratio valid if below                            | 2.0                                                                                                                                               |  |  |

# Analysis log

| image name                               | creation<br>date       | saturation | sampling<br>density | status   |
|------------------------------------------|------------------------|------------|---------------------|----------|
| GreenBeads63x.lif -<br>Series004_bead395 | 2024-10-26<br>09:47:32 | none       | correct             | analysed |

#### Formulas used:

Lateral  $(res_{x,y}^o)$  and axial  $(res_z^o)$  theoretical resolution values used for widefield microscopes are calculated as defined in Wilhelm, S. Confocal Laser Scanning Microscopy, 2011:

$$res_{x,y}^o = \frac{0.51*\lambda_{em}}{NA}$$
  $res_z^o = \frac{1,77n*\lambda_{em}}{NA^2}$ 

NA: numerical aperture,  $\lambda_{em}$ : emission wavelength, n: refractive index of the lens immersion & mounting media.

Axis profiles are fitted using ImageJ Gaussian Curve Fitter and the following formula  $y = a + (b - a) * e^{\frac{-(x-c)^2}{2d^2}}$  (Gaussian fitting).

Measured lateral and axial resolution (Full Width at Half Maximum, FWHM) values are derived using FWHM =  $2d\sqrt{2ln(2)}$ 

Compliance with the Shannon-Nyquist criterion uses the following formulas for Shannon-Nyquist distances calculation:

$$\alpha = \arcsin(\frac{NA}{n})$$

$$\Delta_{x,y} = \frac{\lambda_{em}}{4.NA} \qquad \Delta_z = \frac{\lambda_{em}}{2.n. (1-\cos{(\alpha)})}$$