



63X_1_4_525 - Multiple Bead Image Summary

Microscope info:

Image		Image6				
image's creation	date	2024-10-17 10:22:20				
	method used	from file creation date				
Actual image depth		16				
Microscope type		WideField				
Objective	NA	1.4				
	im. refractive index	1.518				
Channel(s)		Wavelengths		Saturation	sampling (X,Y,Z)	
		Ex. (nm)	Em. (nm)		Nyquist (μm)	Found (μm)
Channel 0			525.0	none	0.094x0.094x0.282	0.063x0.063x0.06
						0.7, 0.7, 0.2

Warnings:

(No saturated pixels detected).

(All channels sampled following Shannon-Nyquist criterion).

(A subresolution bead is used for all channels).


Analysis parameters

Tool & Operator	Tool	Batch PSF Profiler
	Versions	MetroloJ_QC v1.3.1.1, ImageJ v2.14.0/1.54f, Java v1.8.0_322, OS Mac OS X
	Operator & date	SO, October 25, 2024 2:35 PM
data	result folder	/Users/oggsc/Documents/OM/ImageAnalysis/QC/Elyra/PSFs/20241014/63X_1_4/525/
	Type of saved data	.pdf, .jpg, .xls
	Input data bit depth	16
Dimension order		XY-(C)Z
Discard saturated samples		false
Beads	Bead detection threshold	Legacy
	Center detection method	Legacy Maximum Intensity
	Discard bead if more than one particle are thresholded	true
	Background annulus thickness in μm	0.5
	Background annulus distance to bead edges in μm	0.5
	Multiple beads in image	true
	Bead identification method	Using Find Maxima (prominence of 1000.0)
	Bead size (μm)	0.1
	Bead crop Factor	5.0
	Cropped ROI size in μm	2.31x2.31 (using bead size & background annulus parameters)
Square Root PSF Image displayed		true
Tolerance	Applied in this report	true
	X & Y FWHM ratios valid if below	1.5
	Z FWHM ratio valid if below	2.0
Measurement rejected	Outliers	true (using IQR)
	R2 ratio below	0.95

Analysis log

image name	creation date	sampling density	identified raw beads	valid beads	saturation	status
Image 6	2024-10-17 10:22:20	correct	12	10	none	valid beads found
				bead0	none	analysed
				bead1	none	analysed
				bead2	none	analysed
				bead3	none	analysed
				bead4	none	analysed
				bead5	none	analysed
				bead6	none	analysed
				bead7	none	analysed
				bead8	none	analysed
				bead9	none	analysed

Identified beads

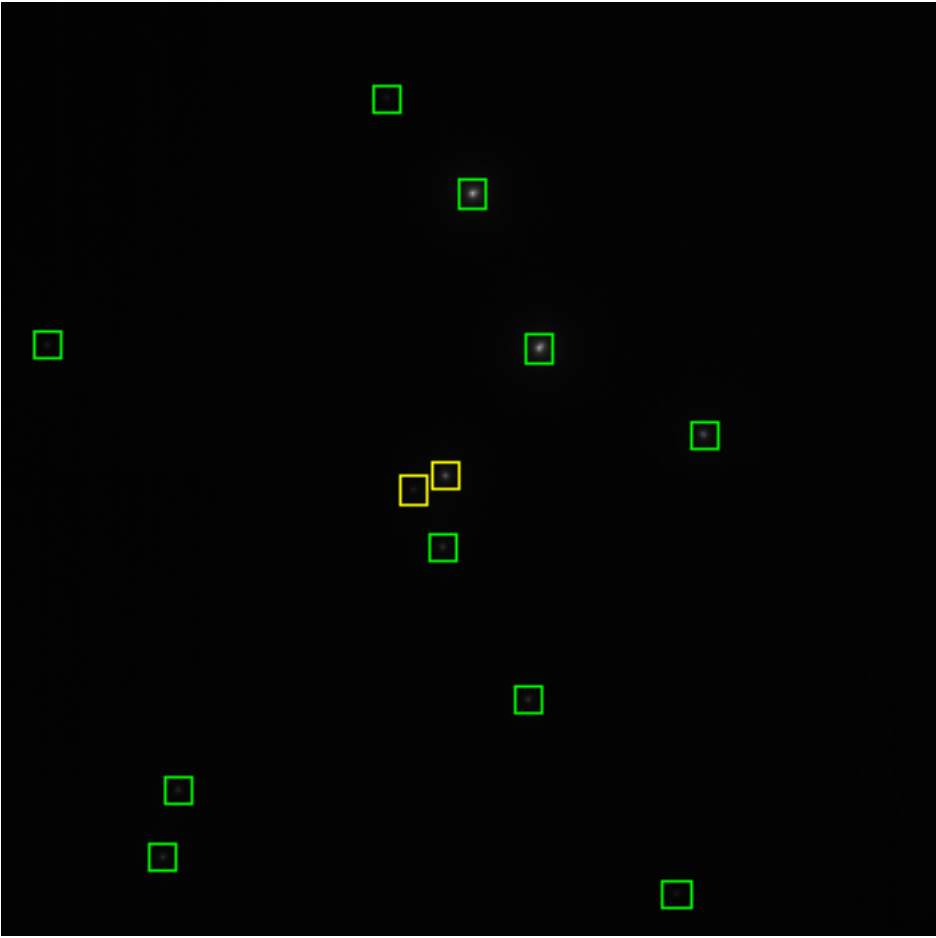

green: valid bead, yellow: too close to another bead, magenta: too close to stack's top or bottom, cyan: too close to the image's edges.