

## **GreenBeadsIMG4R (SHORT)**

### Microscope info:

| Image  |                      | Image6_bead7 |             |            |                       |                      |                     |  |
|--|----------------------|--------------|-------------|------------|-----------------------|----------------------|---------------------|--|
| image's date 2024-10-17 10:22:20             |                      |              |             |            |                       |                      |                     |  |
| creation method used from file creation date |                      |              |             |            |                       |                      |                     |  |
| Actual image depth                           |                      | 16           |             |            |                       |                      |                     |  |
| Microscope type                              |                      | WideField    |             |            |                       |                      |                     |  |
|  | NA                   | 1.4          | 1.4         |            |                       |                      |                     |  |
| Objective                                    | im. refractive index | 1.518        |             |            |                       |                      |                     |  |
|  |                      | Wavel        | ngths       |            | sampling (X,Y,Z)      |                      |                     |  |
| Channel(s)                                   |                      | Ex.<br>(nm)  | Em.<br>(nm) | Saturation | Nyquist (µm)          | Found (µm)           | Nyquist/found ratio |  |
| Channel 0                                    |                      |              | 440.0       | none       | 0.079x0.079x0.<br>236 | 0.063x0.063x0.<br>06 | 0.8, 0.8, 0.3       |  |
| Bead original coordinates(in pixels)         |                      | 63.0, 469.0  |             |            |                       |                      |                     |  |

### Warnings:

(No saturated pixels detected). (All channels sampled following Shannon-Nyquist criterion). (A subresolution bead is used for all channels).

### Resolution table:

| Channel                 | Sig/Backgn<br>d ratio | Dimension | Measured<br>FWHM<br>(µm) | theory (µm) | Fit<br>Goodness | Mes./theory ratio |
|-------------------------|-----------------------|-----------|--------------------------|-------------|-----------------|-------------------|
|                         |                       | Х         | 0.228                    | 0.16        | 1.0             | 1.42              |
| Channel 0 (em. 440.0nm) | 8.1                   | Υ         | 0.235                    | 0.16        | 1.0             | 1.47              |
| 440.01111)              |                       | Z         | 0.686                    | 0.603       | 0.99            | 1.14              |


Green: within specifications, red: outside specifications (ie. XY ratios above 1.5 or Z ratio above 2.0)

#### Lateral asymmetry ratios:

| Channel                 | Ratio |
|-------------------------|-------|
| Channel 0 (em. 440.0nm) | 0.97  |

## Detailed channel detection info:

## Channel #0

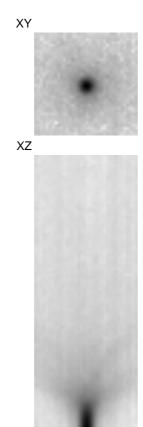



| Channel 0 (em. 440.0nm) |      |           |       |              |  |
|-------------------------|------|-----------|-------|--------------|--|
| Sig./Backgnd ratio      | LAR  | Dimension | FWHM  | Fit goodness |  |
|                         |      | X         | 0.228 | 1.0          |  |
| 8.1                     | 0.97 | Υ         | 0.235 | 1.0          |  |
|                         |      | Z         | 0.686 | 0.99         |  |

### Analysis parameters

|                                 | Tool   | PSF Profiler (batch)   |  |  |
|---------------------------------|--|--|--|--|
| Tool &<br>Operator              | Versions   | MetroloJ_QC v1.3.0, ImageJ v2.14.0/1.54f, Java v1.8.0_322, OS Mac OS X                 |  |  |
|                                 | Operator & date  | aaa, October 20, 2024 9:09 AM  |  |  |
| data                            | result folder  | /Users/bumozaza/Desktop/Zeiss<br>WFM/green/Processed/GreenBeadsIMG4R/Image<br>6/bead7/ |  |  |
| data                            | Type of saved data                                     | .pdf, .jpg, .xls   |  |  |
|                                 | Input data bit depth                                   | 16   |  |  |
| Dim                             | ension order   | XY-(C)Z  |  |  |
| Discard s                       | aturated samples                                       | true   |  |  |
|                                 | Bead detection threshold                               | Legacy   |  |  |
|                                 | Center detection method                                | Centroid   |  |  |
|                                 | Discard bead if more than one particle are thresholded | true   |  |  |
|                                 | Background annulus thickness in µm                     | 0.5  |  |  |
| Beads                           | Background annulus<br>distance to bead edges<br>in µm  | 0.5  |  |  |
| Deads                           | Multiple beads in image                                | true   |  |  |
|                                 | Bead identification method                             | Using Find Maxima (prominence of 1000.0)   |  |  |
|                                 | Bead size (µm)   | 0.1  |  |  |
|                                 | Bead crop Factor                                       | 10.0   |  |  |
|                                 | Cropped ROI size in µm                                 | 2.31x2.31 (using bead size & background annulus parameters)                            |  |  |
|                                 | Bead rejection distance to top/bottom                  | 2.0 μm   |  |  |
| Square Root PSF Image displayed |  | true   |  |  |
|                                 | Applied in this report                                 | true   |  |  |
| Tolerance                       | X & Y FWHM ratios valid if below                       | 1.5  |  |  |
|                                 | Z FWHM ratio valid if below                            | 2.0  |  |  |

# Analysis log

| image name    | creation<br>date       | saturation | sampling density | status   |
|---------------|------------------------|------------|------------------|----------|
| Image 6_bead7 | 2024-10-17<br>10:22:20 | none       | correct          | analysed |

#### Formulas used:

Lateral  $(res_{x,y}^o)$  and axial  $(res_z^o)$  theoretical resolution values used for widefield microscopes are calculated as defined in Wilhelm, S. Confocal Laser Scanning Microscopy, 2011:

$$res_{x,y}^o = \frac{0.51*\lambda_{em}}{NA}$$
  $res_z^o = \frac{1,77n*\lambda_{em}}{NA^2}$ 

NA: numerical aperture,  $\lambda_{em}$ : emission wavelength, n: refractive index of the lens immersion & mounting media.

Axis profiles are fitted using ImageJ Gaussian Curve Fitter and the following formula  $y = a + (b - a) * e^{\frac{-(x-c)^2}{2d^2}}$  (Gaussian fitting).

Measured lateral and axial resolution (Full Width at Half Maximum, FWHM) values are derived using FWHM =  $2d\sqrt{2ln(2)}$ 

Compliance with the Shannon-Nyquist criterion uses the following formulas for Shannon-Nyquist distances calculation:

$$\alpha = \arcsin(\frac{NA}{n})$$

$$\Delta_{x,y} = \frac{\lambda_{em}}{4.NA} \qquad \Delta_z = \frac{\lambda_{em}}{2.n. (1-\cos{(\alpha)})}$$