



## OrangeBeadsReport2

### Microscope info:

Image		Image4_bead4				
image's creation	date	2024-10-17 10:22:29				
	method used	from file creation date				
Actual image depth		16				
Microscope type		WideField				
Objective	NA	1.4				
	im. refractive index	1.518				
Channel(s)		Wavelengths		Saturation	sampling (X,Y,Z)	
		Ex. (nm)	Em. (nm)		Nyquist (μm)	Nyquist/found ratio
Channel 0			590.0	none	0.105x0.105x0.317	0.063x0.063x0.06
Bead original coordinates(in pixels)		480.0, 536.0				

### Warnings:

(No saturated pixels detected). (All channels sampled following Shannon-Nyquist criterion). (A subresolution bead is used for all channels).

### Resolution table:

Channel	Sig/Backgnd ratio	Dimension	Measured FWHM (μm)	theory (μm)	Fit Goodness	Mes./theory ratio
Channel 0 (em. 590.0nm)	5.6	X	0.127	0.215	-0.04	0.59
		Y	2.009	0.215	0.08	9.35
		Z	Center detection failed			

Green: within specifications, red: outside specifications (ie. XY ratios above 1.5 or Z ratio above 2.0)


### Lateral asymmetry ratios:

Channel	Ratio
Channel 0 (em. 590.0nm)	0.06

Detailed channel detection info:

Channel #0

XY

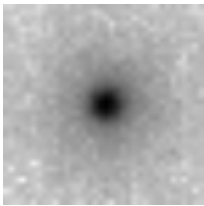

YZ




XZ




Channel 0 (em. 590.0nm)				
Sig./Backgnd ratio	LAR	Dimension	FWHM	Fit goodness
5.6	0.06	X	0.127	-0.04
		Y	2.009	0.08
		Z	Center detection failed	

X profile & fitting parameters:  
Fit equation  $I(x) = a + (b-a) \cdot \exp(-(x-c)/(2 \cdot d))$   
Sum of residuals squared: 588.55462  
Standard deviation: 4.10071  
R^2: 0.07569  
Parameters:  
a = 103.06483  
b = 108.94924  
c = 0.72802  
d = 0.053760

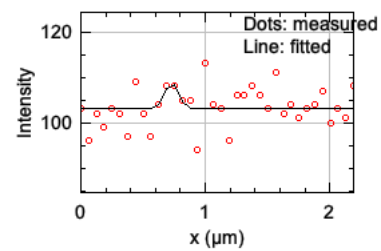

Y profile & fitting parameters:  
Fit equation  $I(y) = a + (b-a) \cdot \exp(-(y-c)/(2 \cdot d))$   
Sum of residuals squared: 759.97502  
Standard deviation: 4.65978  
R^2: 0.18126  
Parameters:  
a = 102.79863  
b = 1.34077E154  
c = -22.64467  
d = 0.85307

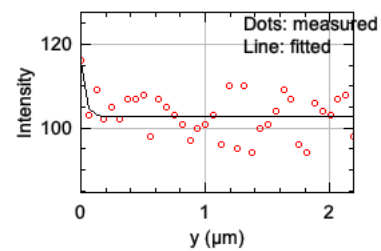

Analysis parameters

Tool & Operator	Tool	PSF Profiler (batch)
	Versions	MetroloJ_QC v1.3.0, ImageJ v2.14.0/1.54f, Java v1.8.0_322, OS Mac OS X
	Operator & date	aaa, October 20, 2024 10:52 AM
data	result folder	/Users/bumozaza/Desktop/Zeiss WFM/orange/Processed/OrangeBeadsReport2/Image 4/bead4/
	Type of saved data	.pdf, .jpg, .xls
	Input data bit depth	16
Dimension order		XY-(C)Z
Discard saturated samples		true
Beads	Bead detection threshold	Legacy
	Center detection method	Centroid
	Discard bead if more than one particle are thresholded	true
	Background annulus thickness in $\mu\text{m}$	0.5
	Background annulus distance to bead edges in $\mu\text{m}$	0.5
	Multiple beads in image	true
	Bead identification method	Using Find Maxima (prominence of 1000.0)
	Bead size ( $\mu\text{m}$ )	0.1
	Bead crop Factor	10.0
	Cropped ROI size in $\mu\text{m}$	2.31x2.31 (using bead size & background annulus parameters)
	Bead rejection distance to top/bottom	2.0 $\mu\text{m}$
Square Root PSF Image displayed		true
Tolerance	Applied in this report	true
	X & Y FWHM ratios valid if below	1.5
	Z FWHM ratio valid if below	2.0

Analysis log

image name	creation date	saturation	sampling density	status
Image 4_bead4	2024-10-17 10:22:29	none	correct	analysed

### Formulas used:

Lateral ( $res_{x,y}^o$ ) and axial ( $res_z^o$ ) theoretical resolution values used for widefield microscopes are calculated as defined in Wilhelm, S. Confocal Laser Scanning Microscopy, 2011:

$$res_{x,y}^o = \frac{0.51 \cdot \lambda_{em}}{NA} \quad res_z^o = \frac{1.77 n \cdot \lambda_{em}}{NA^2}$$

NA: numerical aperture,  $\lambda_{em}$ : emission wavelength, n: refractive index of the lens immersion & mounting media.

Axis profiles are fitted using ImageJ Gaussian Curve Fitter and the following formula  $y = a + (b - a) * e^{\frac{-(x-c)^2}{2d^2}}$  (Gaussian fitting).

Measured lateral and axial resolution (Full Width at Half Maximum, FWHM) values are derived using  $FWHM = 2d\sqrt{2\ln(2)}$

Compliance with the Shannon-Nyquist criterion uses the following formulas for Shannon-Nyquist distances calculation:

$$\alpha = \arcsin\left(\frac{NA}{n}\right)$$

$$\Delta_{x,y} = \frac{\lambda_{em}}{4 \cdot NA} \quad \Delta_z = \frac{\lambda_{em}}{2 \cdot n \cdot (1 - \cos(\alpha))}$$